Bright™

027955,
027956

ELECTRIC TOOTHBRUSH

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Jula AB reserves the right to make changes to the product. Jula AB claims copyright on this documentation. It is not allowed to modify or alter this documentation in any way and the manual shall be printed and used as it is in relation to the product. For the latest version of operating instructions, refer to the Jula website.

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**JULA AB
BOX 363, 532 24 SKARA, SWEDEN**

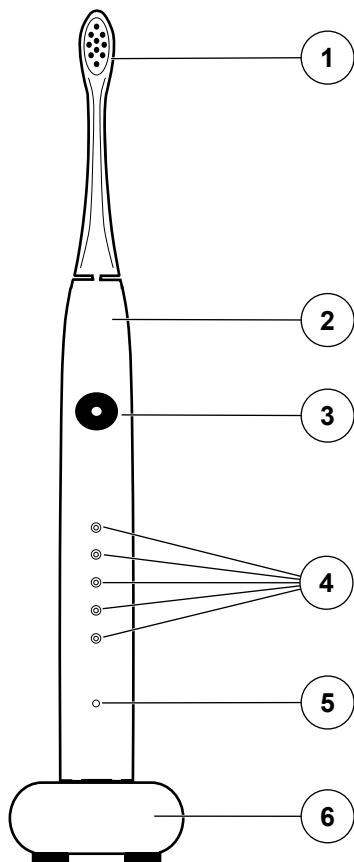
Table of contents

1	Introduction	5
1.1	Product description	5
1.2	Symbols	5
1.3	Product overview	5
2	Safety	6
2.1	Safety definitions	6
2.2	Safety instructions for operation	6
3	Operation	7
3.1	To charge the product	7
3.2	To operate the product	7
4	Storage	8
5	Maintenance	9
5.1	To clean the product	9
6	Disposal	9
6.1	To discard the product	9
7	Technical data	10

1 Introduction

1.1 Product description

The product is an electric toothbrush with 5 different modes. The product is not recommended to be used by children under 3 years of age. Only for domestic use.

1.2 Symbols

	Read the operating instructions carefully and make sure that you understand the instructions before you use the product. Save the instructions for future reference.
	The product has a detachable power supply unit. Only use the detachable supply unit provided with the product.
	Double insulation (safety class II), does not need to be earthed.
	This product complies with applicable EU directives and regulations.
	Recycle as electrical waste.

1.3 Product overview

(Figure 1)

1. Brush head
2. Handle
3. On/off button
4. Mode indicators
5. Charge indicator
6. Charger

2 Safety

2.1 Safety definitions

Warning! If you do not obey these instructions, there is a risk of death or injury.

Caution! If you do not obey these instructions, there is a risk of damage to the product, other materials or the adjacent area.

Note! Information that is necessary in a given situation.

2.2 Safety instructions for operation

Warning! If you do not obey these instructions, there is a risk of death or injury.

- Read the warning instructions that follow before you use the product.
- Do not modify the product.
- The product must only be used for its intended purpose.
- Do a check for damage on the product, the power cable and the accessories regularly. Do not use the product if there is any sign of damage. Refer to an authorized service center for reparation.
- The product can be used by children over 3 years of age and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Do not let children use the product without adult supervision.
- Cleaning and maintenance shall not be done by children.
- Do not let children play with the product.
- Do not use accessories that are not adapted for the product.
- Do not put the charger or the power plug into water or any other liquid, or where it can fall or be pulled into a tub or sink.

- If the charger falls into the water or any other liquid, pull out the power plug immediately. Do not reach for the charger when it is in water until the power plug has been pulled out.
- The battery of the product is non-replaceable.
- Do not attempt to disassemble the product except when the battery is to be disposed.
- Do not pull the power cable to unplug the power plug.
- Do not touch the power plug with wet hands, this can cause an electric shock.
- If you are under treatment for any oral care condition, consult your dental professional before you use the product.
- Follow the instructions for maintenance to prevent brush head breakage that can cause choking hazard or damage the teeth. Refer to "5 Maintenance" on page 9.
- Do not use the product without a brush head.

3 Operation

3.1 To charge the product

Warning! Only a charger intended for the product must be used.

- 1 Attach the charger into a power outlet.
- 2 Put the handle on the charger.

3.2 To operate the product

Note! Your gums can bleed slightly the first few times you use the product. In general, this will stop after a few days. If it persists after 2 weeks, refer to a dentist or dental hygienist.

Note! You can experience a ticking or tingling sensation the first few times you use the product, this will subside after a period of time. To alleviate the sensation, use the sensitive mode.

- 1 Put the brush head on the handle.
- 2 Wet the brush head and apply toothpaste on it.

- 3 Put the brush head on the outside surface of one of your teeth and then start the product.

Note! The bristles of the brush head should be placed against the teeth at an angle towards the gum line.

- 4 Brush the surface of the teeth for a few seconds and then move the brush head to the next teeth.
- 5 Brush the inside and outside surfaces of your teeth throughout the mouth for the recommended 2 minutes and then brush the chewing surfaces of your teeth.
- 6 Before you remove the brush head from the mouth, push the on/off button to stop the product.

3.2.1 Brushing modes

Mode	Description
Clean	Regular tooth brushing.
Sensitive	Brushing for sensitive areas.
White	Brushing to remove surface stains.
Polish	Brushing to increase the shine of the teeth.
Massage	Massage brushing that is gentle to the gums.

3.2.2 Brushing intervals

The product does a short pause at 30 second intervals as a help to brush the 4 quadrants of the mouth equally. A long pause indicates that the recommended brushing time of 2 minutes has been reached.

4 Storage

The handle can be stored on the charger. The charger can be kept plugged into the power outlet to maintain full power if the product is used everyday.

Note! For environmental reasons it is recommended to disconnect the charger until the battery is discharged.

5 Maintenance

Warning! The product has no parts that can be serviced by the user. Refer to an authorized service center for repairation.

- Make sure that the brush head fits well on the handle before each use. If the brush head does not fit well, replace the brush head before use.
- If you dropped the product, replace the brush head before use even if there is no visible damage to it.
- Replace the brush head every 3 months or when the brush head shows signs of wear.
- To maintain the maximum capacity of the battery, disconnect the charger and let the battery become fully discharged at least every 6 months.

5.1 To clean the product

- Clean the brush head after each use.
- Clean the charger with a damp cloth regularly.
- Regularly rinse the brush head and the handle thoroughly under running water. Make sure to wipe the brush head and the handle dry before you put the handle back on the charger.

6 Disposal

6.1 To discard the product

- Make sure that you follow local regulations when you discard the product. Do not burn the product.

7 Technical data

Specification	Value
Rated voltage	100-240 / 3.7 V
Battery	Li-ion
Movement	38000 /min
Charging time	6 h
Charger	Z57

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Jula AB, Box 363, SE-532 24 SKARA, SWEDEN