MASSAGE GUN

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Jula AB reserves the right to make changes to the product. Jula AB claims copyright on this documentation. It is not allowed to modify or alter this documentation in any way and the manual shall be printed and used as it is in relation to the product. For the latest version of operating instructions, refer to the Jula website.

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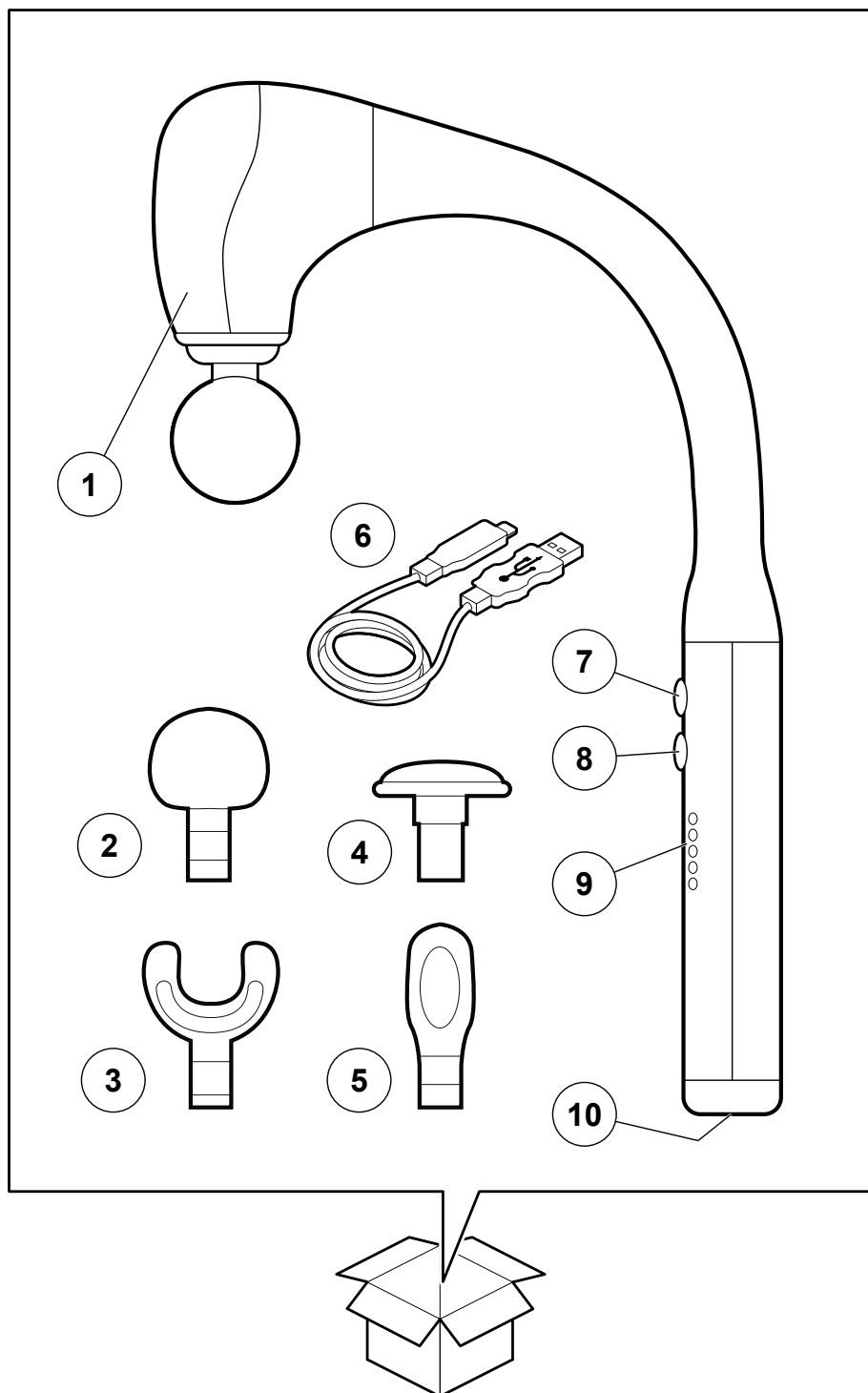
1 Introduction

1.1 The product

The product is a U-shaped massage gun with 5 different massage modes. 4 different massage heads are included.

1.2 Dedicated use

The product is used to massage sore and stiff body parts. The shape of the product makes it easier to massage the back and shoulder blades. Intended for household use.

1.3 Symbols

	Read the operating instructions carefully and make sure that you understand the instructions before you use the product. Save the instructions for future reference.
	This product complies with applicable EU directives and regulations.
	Recycle as electrical waste.

1.4 Product overview

1. Massage gun
2. Round head
3. Fork head
4. Flat head
5. Cylindrical head
6. USB-C cable
7. Power button
8. Mode button
9. Mode/ battery display indicator
10. Charging port

2 Safety

2.1 Safety definitions

⚠ Warning! If you do not obey these instructions, there is a risk of death or injury.

⚠ Caution! If you do not obey these instructions, there is a risk of damage to the product, other materials or the adjacent area.

Note! Information that is necessary in a given situation.

2.2 Safety instructions for operation

⚠ Warning! Stop using the product immediately if you experience pain or discomfort.

- Read and obey the warning instructions before operation.
- Make sure that these persons do not operate the product independently:
 - Children less than 8 years.
 - Persons who do not know the product sufficiently.
 - Persons with disabilities.
- Children must not clean or do maintenance on the product independently.
- Do not use the product if it is damaged.
- Do not make changes to the product.
- Use the product only for its specified function.
- Do not cover the product or the air flow during use, the product can become hot.
- Keep the product away from water and do not touch it with wet hands.
- Only use the correct adapter to charge the product, refer to "8 Technical data" on page 6.
- Keep the product monitored when it charges.
- Only use the included accessories.

- Do not operate the product when it charges.
- Use only on soft body tissues. Do not use the product if it cause pain or discomfort.
- Do not use the product on hard areas of the body or near sensitive areas. For example on the head, bones, shins, on the top of the feet and hands, eyes, breasts and genitals.
- To use the product safely:
 - Do not push too hard. Too much pressure can cause bruises.
 - Look at the massage area often.
 - If you feel pain or discomfort, stop immediately.
- Only use the product on dry, clean skin.
- Keep your fingers, hair and other parts of the body away from the axle and the back of the massage head because of the pinch risk.
- Do not use the product when you walk or operate a vehicle.
- Do not use the product in a way that can damage bones.
- Do not use the product in moist conditions, or at temperatures over 40°C.
- Be careful when you use the product if you are pregnant, have a peacemaker, have high blood pressure, are injured, are under medical treatment or are sick.
- If there is continuous pain after usage, speak to a doctor.

3 Operation

3.1 To do before you operate the product

- Charge the product for 3 hours before you use it for the first time.
- Remove jewellery and other accessories that can have an unwanted effect on the massage.

3.2 Massage gun-head overview

1. **Round head:** For large muscle groups.
2. **Fork head:** For soft tissues around the neck, spine and calves.
3. **Flat head:** For all body parts.
4. **Cylindrical head:** For joints and to hit deep tissues.

3.5 To change massage heads

- 1 Pull out the attached gun-head.
- 2 Attach the new gun-head.

3.6 To operate the product

Note! The product will automatically stop after 15 min without operation.

Note! The product has an overheating protection.

The product has 5 different massage modes.

- 1 Select the applicable massage head.
- 2 Push the power button for 2 seconds to start the product.
- 3 Push the mode button for 1 second to change between the different modes.
- 4 Push the power button for 2 seconds to stop the product.

3.6.1 The battery indicator

Note! When you charge the product the mode/ battery indicator will flash gradually and you will not be able to use the product.

- **All the mode/ battery indicator stay on:** the battery is fully charged.
- **All the mode/ battery indicator flashes for 10 seconds:** the battery is low.

4 Storage

- Keep the product in a dry area, protected from direct sunlight and wet conditions, where children and pets do not have access to it.

5 Maintenance

 **Warning!** Make sure the product is stopped before you do any maintenance.

 **Warning!** Do not put or wash the product in water.

 **Caution!** Do not use strong corrosive agents for example: petrol, or thinner to clean the product.

- Clean the product with a damp cloth.
- If you do not use the product for a long time, charge it every 2 months for the preservation of the battery capacity.

6 Troubleshooting

Problem	Possible cause	Task
The product does not start.	The battery level is low.	Make sure the product is connected to a serviceable power outlet. Charge the product for 3 hours.
The product stops during use.	The overheat protection activated.	Let the product become cool before you start it again.
	The battery level is low.	Charge the product for at least 3 hours.

7 Discard

- Make sure that you follow local regulations when you discard the product. Do not burn the product.

8 Technical data

Specification	Value
Voltage	230 V ~50 Hz
Output	5 V
Battery capacity	1500 mAh
Power consumption in off mode (in watts rounded to the first decimal place)	0.13 W
Adapter	RSS1002-100050-W2E
Number of speed modes	5
Battery charging time	3 h
Cord length	50 cm
Dimensions (LxWxD)	34.5 x 20 x 6 cm
Weight	481 g

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