MASSAGE GUN HOT AND COLD

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Jula AB reserves the right to make changes to the product. Jula AB claims copyright on this documentation. It is not allowed to modify or alter this documentation in any way and the manual shall be printed and used as it is in relation to the product. For the latest version of operating instructions, refer to the Jula website.

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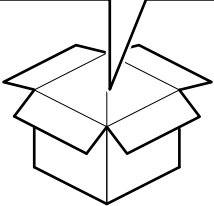
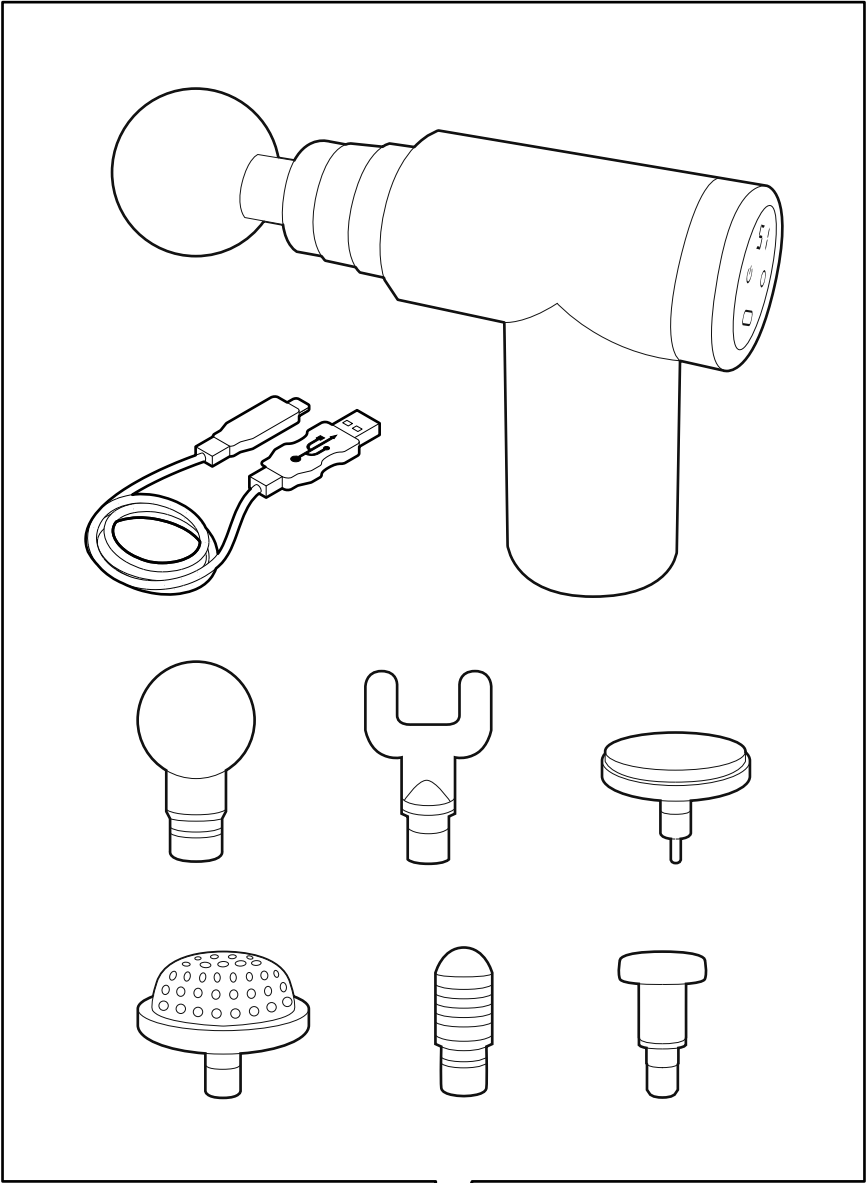
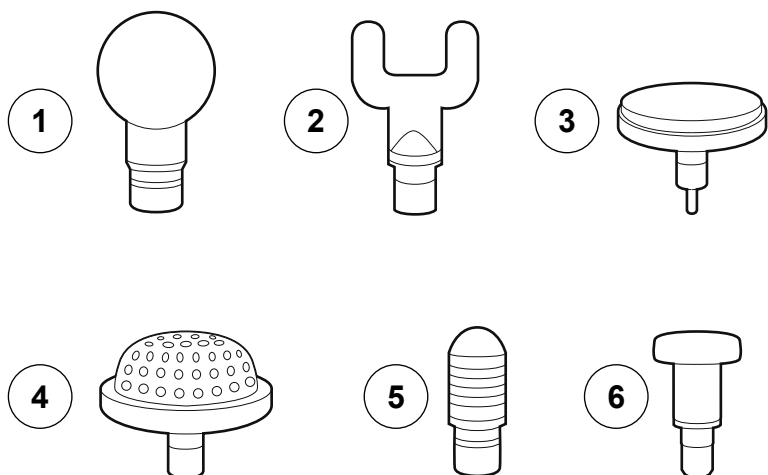
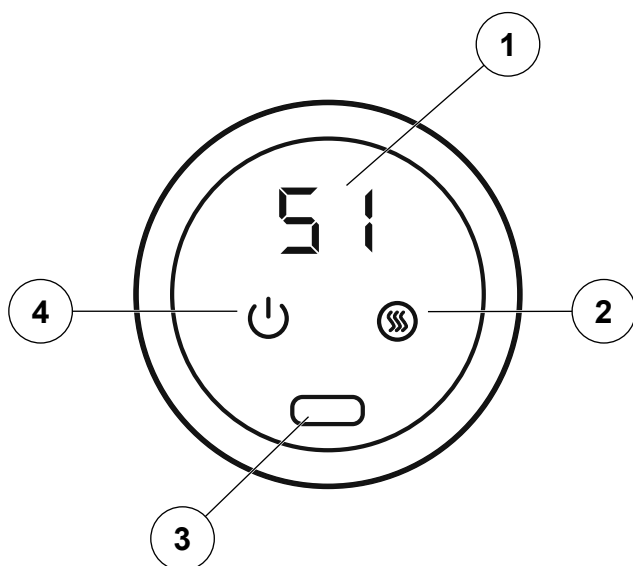
1 Introduction

1.1 The product

The product is a massage gun with heat and cool function. There are 6 different massage heads included.

1.2 Dedicated use

The product is used to massage sore and stiff body parts. Intended for household use.

1.3 Symbols

	Read the operating instructions carefully and make sure that you understand the instructions before you use the product. Save the instructions for future reference.
	This product complies with applicable EU directives and regulations.
	Recycle as electrical waste.

2 Safety

2.1 Safety definitions

Warning! If you do not obey these instructions, there is a risk of death or injury.

Caution! If you do not obey these instructions, there is a risk of damage to the product, other materials or the adjacent area.

Note! Information that is necessary in a given situation.

2.2 Safety instructions for operation

Warning! Stop using the product immediately if you experience pain or discomfort.

- Read and obey the warning instructions before operation.
- Do not use the product if it is damaged.
- Do not make changes to the product.
- Use the product only for its specified function.
- Make sure that these persons do not operate the product independently:
 - Children less than 8 years.
 - Persons who do not know the product sufficiently.
 - Persons with disabilities.
- Children must not clean or do maintenance on the product independently.
- Do not cover the product or the air flow during use, the product can become hot.
- Keep the product away from water and do not touch it with wet hands.
- Only use the correct adapter to charge the product, refer to "8 Technical data" on page 7.
- Keep the product monitored when it charges.
- Only use the included accessories.
- Do not operate the product when it charges.
- Use only on soft body tissues. Do not use the product if it cause pain or discomfort.
- Do not use the product on hard areas of the body or near sensitive areas. For example on the head, bones, shins, on the top of the feet and hands, eyes, breasts and genitals.

- Only use the product on dry, clean skin.
- To use the product safely:
 - Do not push too hard. Too much pressure can cause bruises.
 - Look at the massage area often.
 - If you feel pain or discomfort, stop immediately.
- Keep your fingers, hair and other parts of the body away from the axle and the back of the massage head because of the pinch risk.
- The product has a heated surface. Be careful if you have sensitive skin.
- Do not use the product in moist conditions, or at temperatures over 40°C.
- Do not use the product when you walk or operate a vehicle.
- Be careful when you use the product if you are pregnant, have a pacemaker, have high blood pressure, are injured, are under medical treatment or are sick.
- Do not use the product in a way that can damage bones.
- If there is continuous pain after usage, speak to a doctor.

3 Operation

3.1 To do before you operate the product

- Charge the product for 3 hours before you use it for the first time.
- Remove jewellery and other accessories that can have an unwanted effect on the massage.

3.2 Massage gun-head overview

Figure 1

1. **Round head:** For large muscle groups.
2. **Fork head:** For soft tissues around the neck, spine and calves.
3. **Hot compress head:** For stiff neck and to relieve muscle pain.
4. **Cold compress head:** For joint swelling and muscle inflammation.
5. **Cylindrical head:** For joints.
6. **Flat head:** For all body parts.

3.7 Display overview

Figure 2

1. Display for battery capacity/ mode
2. Heat mode button
3. Charging port
4. Power/ mode button

3.5 To change massage heads

- 1 Pull out the attached gun-head.
- 2 Attach the new gun-head.

3.6 To operate the product

Note! The product will automatically stop after 15 min without operation.

Note! The product stops automatically when the battery level is 5% or below. Charge the battery before operation.

Note! The product has an overheating protection.

The product has 6 different massage modes.

- 1 Select the applicable massage head.
- 2 Push the power button for 2 seconds to start the product.
- 3 Push the power button repeatedly to change between massage modes.
- 4 Push the power button for 2 seconds to stop the product.

3.6.1 To use the heating mode

- 1 Attach the hot-compress head.

- 2 Push the heating mode button for 2 seconds to start the heat mode.

- 3 Push the power button and select the applicable massage mode.

3.6.2 To use the cold mode

Note! Put the cold-compress head massage accessory in the freezer for 2 hours before use.

- 1 Attach the cold-compress head.
- 2 Push the power button and select the applicable massage mode.

6 Troubleshooting

Problem	Possible cause	Task
The product does not start.	The battery level is low.	Make sure the product is connected to a serviceable power outlet. Charge the product for 3 hours.
The product stops during use.	The overheat protection activated.	Let the product become cool before you start it again.
	The battery level is under 5%.	Charge the product for at least 3 hours.

7 Discard

- Make sure that you follow local regulations when you discard the product. Do not burn the product.

4 Storage

- Keep the product in a dry area, protected from direct sunlight and wet conditions, where children and pets do not have access to it.

5 Maintenance

Warning! Make sure the product is stopped before you do any maintenance.

Warning! Do not put or wash the product in water.

Caution! Do not use strong corrosive agents for example: petrol, or thinner to clean the product.

- Clean the product with a damp cloth.
- If you do not use the product for a long time, charge it every 2 months for the preservation of the battery capacity.

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