MASSAGE GUN

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SAFETY INSTRUCTIONS

Read these instructions carefully before use, and save them for future reference. Failure to follow the instructions below can result in a risk of burn injuries, fire, electric shock and/or personal injury.

- Only use the product on soft body tissue and in such a way that it does not cause pain or discomfort.
- Do not use the product on the head, skeleton bones, or hard parts of the body.
- Moderate or excessive pressure can result in bruising. Check the massaged area at regular intervals. Stop using the product immediately if you experience pain or discomfort.
- Only use the product on dry, clean skin. Press gently and move the product over the skin.
- Keep your fingers, hair and other parts of the body away from the axle and the back of the massage head – pinch risk.
- Do not use the product where skeleton bones are covered by thin layers of muscles, such as on the shins, the top of the feet or the back of the hands.
- Do not use the product on or near sensitive areas, such as the eyes, teeth, breasts, genitals, artificial limbs, etc.
- Do not use the product in such a way as to risk fracture.
- In the event of persistent pain, stop using the product immediately and seek medical attention.
- Always check the product, the power cord and charger before use. If the power cord or charger get damaged they must be replaced by an authorised service centre.
- This product can be used by children from 8 years and upwards and by persons with physical, sensorial or mental disabilities, or persons who lack experience or knowledge, if they are supervised or receive instructions concerning the

safe use of the product and understand the risks involved with its use. Do not allow children to play with the product. Do not allow children to clean or maintain the product without supervision.

HOW TO USE

- Do not press the product too hard against soft tissue – risk of personal injury. Excessive pressure activates the safety function and switches off the product. If this happens, wait 3 seconds and then start the product again.
- Never poke any objects in the ventilation openings on the product.
- Do not place or store the product where it could drop, or be pulled down into a bath or washbasin.
- Never immerse the product in water or any other liquid.
- Never use the product under a quilt, blanket or cushion. The product can overheat – risk of fire, electric shock and/or personal injury.
- Never modify the product in any way.
- Do not start the product when charging is in progress.
- Do not expose the product to direct sunlight.
- Do not use the product for more than 15 minutes at a time. Allow the product to cool for 30 minutes after using it for 15 minutes before using it again.

WARNING!

- **Do not touch the product, power switch, charger or plug with wet hands.**
- **Switch off the product after use.**
- **Never use the product if you are tired or under the influence of drugs, alcohol or medication.**
- **Do not use the product in wet or explosive environments, for example in the vicinity of flammable liquids, gas or dust. The product produces sparks, which can ignite dust and fumes.**

- **Do not use the product in damp conditions, or at temperatures over 40°C.**

IMPORTANT:

- **The product is only intended for household use.**
- **Always consult a doctor before using the product if you are pregnant, wear a pacemaker, or have issues concerning your health and the use of this product.**
- **Never cover the product when in use.**
- **Do not use the product for more than 15 minutes at a time.**
- **Do not use the product immediately before going to bed.**
- **Prolonged use can cause the product to overheat and reduce its life span. If the product gets too hot, switch off the product, pull out the plug, and allow it to cool before using it again.**
- **Keep the product under careful supervision when it is used by, or near children or persons with functional disorders.**

- **Only use accessories and attachments recommended by the manufacturer.**

BATTERIES

- Only use the supplied charger, otherwise there is a risk of electric shock, fire and/or material damage.
- Do not expose batteries to water or any other liquid.
- Do not charge batteries in high temperatures.
- Do not use batteries near naked flames or other sources of heat. Move the battery to a cooler place immediately if it leaks or smells.
- Never short-circuit the battery terminals.
- When batteries are not in use keep them away from metal objects such as paperclips, coins, keys, nails, screws and other small metal objects that can short circuit the terminals and cause a fire and/or explosion.
- Do not puncture batteries, tread on them or subject them to knock and blows.

- Never attempt to open or dismantle batteries.
- The battery must be removed before the product is handed in for recycling.
- Switch off the product before removing the battery.
- Recycle used batteries in accordance with local regulations.

CHARGING


Make sure to do the charging in a suitable place and to use the right charger. Never connect an unknown charger, or a charger, charging cord or battery that can be suspected to be defective. Otherwise there is a risk of overheating, and in the worst case fire. Do not leave the battery unattended when charging.

See website of Swedish National Electrical Safety Board for further information.

WARNING!

Only charge the battery with the supplied charger (RSS1002-126126-W2E).

SYMBOLS

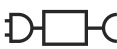	Read the instructions.
	Only use the supplied removable plug.
	Approved in accordance with the relevant directives.
	Recycle discarded product in accordance with local regulations.

TECHNICAL DATA

Rated voltage	12 VDC
Power consumption in off mode (in watts rounded to the first decimal place)	0.09 W
Stroke rate	0-3200/min
Battery technique	Li-ion
Battery capacity	2000 mAh
Battery life	6-8 h
Charging time	2-5 h
Noise level	35-45 dB
Weight	0.9 kg

DESCRIPTION

Portable, cordless massager for tender muscles and joints, or for preventive massage. With six replaceable massage heads and 30 different speeds and modes. Also has digital display, battery level indicator, non-slip handle and radiator. Charger included. Parts:

1. *Massage head*
2. *Opening for cool air*
3. *Display*
4. *Handle with non-slip rubber grip*
5. *Charging connector*
6. *Power switch*

FIG. 1

USE

CHARGING

Charge the battery for 4 hours before using for the first time. When the battery is fully charged the status light goes from red to green.

- The power switch must be in the off position when charging the battery.
- Make sure that charging area is well ventilated and free from foreign objects, water and any other liquids, and that there are no sources of heat nearby.
- Do not use the product when charging is in progress.
- Battery is inserted in the handle and it is normal for the handle and charger to become warm when charging.
- When the battery is not in use keep it away from metal objects such as paper-clips, coins, keys, nails, screws and other small metal objects that can short circuit the terminals.
- Stop charging and/or using immediately if the battery becomes deformed, abnormally warm, smells, leaks, or otherwise out of order.

HOW TO USE

1. Press one of the six massage heads on the massage head holder and twist to lock in place.
2. Switch on the power switch on the bottom of the handle. The status light goes on to show that the power is switched on.
3. Press the power switch symbol on the display.
4. The battery charge level is shown in the window on the right: 99 for fully charged battery, dropping to 80, 60, 40 and 20 as the battery discharges.
5. The speed is shown in the left window. Press the + and – buttons repeatedly to step through the speeds until the required speed is set.
6. Switch off the product after use.

NOTE:

- **The product must only be used on dry skin.**
- **Do not use skin care products, such as skin lotion or oils.**
- **The massage heads should not be rinsed or washed.**
- **Do not use the product on head, abdomen or joints. Incorrect use can cause discomfort (consult doctor before use).**

DIFFERENT MASSAGE HEADS

Spherical massage head

Recommended for all groups of muscles – arms, hips, legs etc.

FIG. 2

Stroke the massage head slowly along the muscles and let the vibrations relax them. Hold the massage head still at places that ache or are tender.

U-shaped massage head

Recommended for the spine.

FIG. 3

Do not use the U-shaped massage head directly on vertebrae. Move the massage head slowly along the muscles and allow the vibrations to relax tense muscles. Hold the massage head still at places that ache or are tender.

Projectile shaped massage head

Recommended for acupuncture areas and the soles of the feet.

FIG. 4

Stroke the massage head slowly along the muscles and let the vibrations relax them. Hold the massage head still at places that ache or are tender.

Flat circular massage head

Recommended for full body massage – neck, back, tendons etc.

FIG. 5

Stroke the massage head slowly along the muscles and let the vibrations relax them. Hold the massage head still at places that ache or are tender.

Arch shaped massage head

Recommended for shoulders and legs.

FIG. 6

Stroke the massage head slowly along the muscles and let the vibrations relax them. Hold the massage head still at places that ache or are tender.

Flat massage head with pressure points

Recommended for full body massage.

FIG. 7

Stroke the massage head slowly along the muscles and let the vibrations relax them. Hold the massage head still at places that ache or are tender.

MAINTENANCE

- Clean the product with a damp cloth and wipe dry with a dry cloth.
- Do not use benzene, thinner or other strong or corrosive detergents.
- Always check before use that the massage heads are intact and firmly in place. Loose or damaged parts can affect functionality.
- When the product is not in use it should be stored in its case in a clean, dry and well ventilated place. Do not expose the product to moisture, extreme temperatures or direct sunlight.

TROUBLESHOOTING

Follow the troubleshooting instructions below in the event of any problems with this product.

Problem	Possible cause	Remedy
The product switches off when in use.	Battery temperature too high.	Restart the product.
	After 20 minutes of continuous use the overheating protection for the battery switches off the product	Restart the product.
	The massage head has been pressed too hard against the body. The overload protection for the motor then switches off the motor.	Press the + button to restart the motor.

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Jula AB reserves the right to make changes to the product. Jula AB claims copyright on this documentation. It is not allowed to modify or alter this documentation in any way and the manual shall be printed and used as it is in relation to the product. For the latest version of operating instructions, refer to the Jula website.