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Bruksanvisning for TENS-apparat

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SAFETY INSTRUCTIONS

Read the User Instructions carefully before use.

The following symbols are used in these User Instructions to highlight important instructions.

Symbol	Meaning
 DANGER!	Indicates an immediately dangerous situation that can lead to fatal or serious personal injury if the instructions are not followed.
 WARNING!	Indicates a potentially dangerous situation that could lead to personal injury and/or material damage if the instructions are not followed.
 IMPORTANT!	Indicates a potentially dangerous situation that could lead to material damage if the instructions are not followed.

Symbol	Meaning
	Indicates prohibited measures, e.g. dismantling the unit.
	Indicates instructions that must be followed.
	BF type, patent applied for part.

 DANGER!	
Do not use the unit together with other medical electronic equipment.	
 WARNING!	
If you are pregnant, consult a doctor before using the unit.	
<ol style="list-style-type: none"> 1. Do not place the electrodes over the heart, head or mouth, or on the genitals or injured skin. 2. Do not use the unit in high temperatures or in the vicinity of flammable liquids, gases or dust. 3. Do not use the unit in the bath or shower, or when sleeping. 	
Do not dismantle the unit, or attempt to modify or repair it.	

IMPORTANT!

1. Do not move the unit when in use.
2. Keep the unit under careful supervision when it is used by, or near children, or persons with functional disorders.
3. Stop using the unit immediately and seek medical attention if you experience discomfort during use.
4. Remove the batteries from the unit if it is not going to be used for some time.

TECHNICAL DATA

Electrical data	4.5 VDC, 100 mA
Display	Digital LCD display
Operating conditions:	
- Temperature	5 to 40°C
- Relative air humidity	≤ 80%
- Atmospheric pressure	86–106 kPa
Electrical output data:	
- Pulse frequency	1–1000 Hz
- Pulse time	20–400 μs
- Impedance	500 Ω, precision ± 10%
- Treatment time	15 minutes, precision ± 10%
Transport and storage:	
- Ambient temperature	-20 to 55°C
- Relative air humidity	≤ 93%
- Atmospheric pressure	50–106 kPa

DESCRIPTION


The unit is microprocessor controlled and combines traditional Chinese massage with modern technology. There are four different modes to stimulate tissues and promote the blood circulation.

Function

The unit stimulates tissue with the help of weak electric current at low frequency, which expands the blood vessels, increases the blood circulation and facilitates the removal of waste products. At higher current strengths the unit can stimulate nerves and cause muscular contraction, which counteracts fatigue and can contribute to the treatment of peripheral nerve damage and paralysis.

Functions and features

- New design, easy to use.
- 6 automatic and 4 manual modes to suit everyone.
- Increases the blood circulation, counteracts muscular stiffness, inflammation, swelling, pain and fatigue. Effective and convenient.
- Electronic stimulation of muscles and tissues.
- Large, comprehensible display shows active function.
- Select the required area, depending on symptoms.

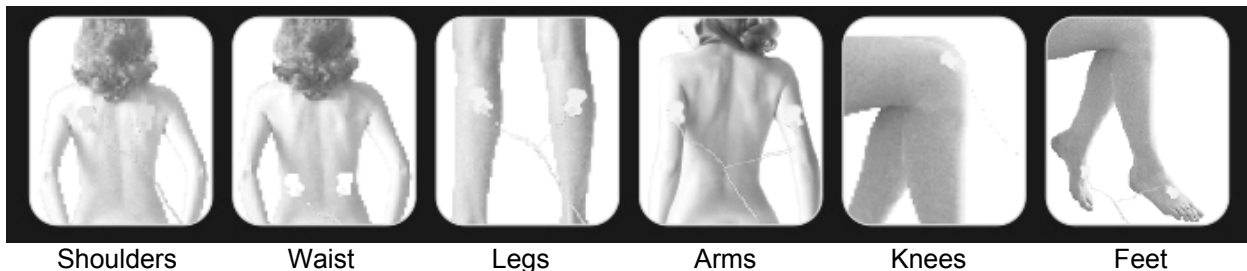Shoulders

Waist


Legs

Arms

Knees

Feet

Parts

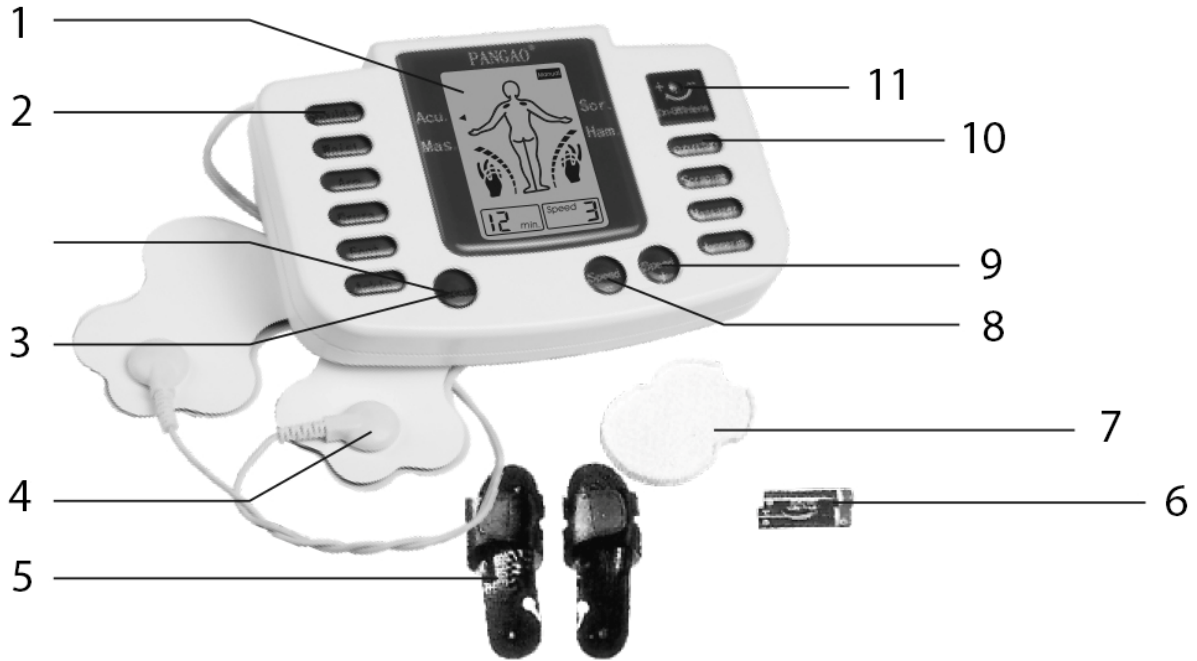- | | |
|------------------------------------|---------------------------------|
| 1. LCD display | 7. Fastener |
| 2. Automatic mode med 6 settings | 8. Reduction of pulse frequency |
| 3. Repeat | 9. Increase of pulse frequency |
| 4. Electrode | 10. Manual mode med 4 settings |
| 5. Foot massage accessory (option) | 11. Switch/intensity setting |
| 6. Battery | |

USE

Preparations

1. Power supply
 - Insert 3 x AA batteries in the battery compartment.
 - Connect the mains adapter to the unit and then plug the adapter into a power point.
2. Connect the electrode leads.

Adjust to the required mode

1. Auto: 6 automatic modes
2. Manual: 4 manual modes
3. Speed+: higher pulse frequency
4. Speed-: lower pulse frequency
5. Repeat: locking of the selected mode

Starting

1. Remove the protective film from the electrodes and attach them to the body at the required place.
2. Start the unit with the power switch (On/+).

Use


1. Press the power switch (On/+). Automatic mode 1 (shoulder massage 15 minutes) is shown on the display.
2. Select the required intensity from 1 to 8. NOTE: Change the intensity gradually.
3. Modes
4. Automatic modes for shoulders, waist, legs, arms, knees and feet. 4 manual treatment modes.
5. Set the required pulse frequency with the Speed+ and Speed- buttons. 1 is the lowest and 8 is the highest.
6. The set mode can be locked with the Repeat button in automatic mode. Press Repeat again to reactivate switching between modes.

NOTE: Switch off the unit after use and store in a dry and safe place.


Positioning the electrodes for common complaints

Shoulder inflammation

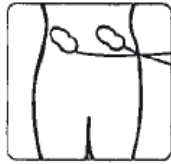Chronic back complaints

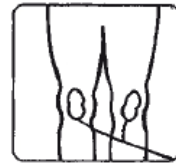Arthritis

Accessory for foot massage

1. Connect the electrode leads to the unit and to the foot massage unit.
2. Put both your feet in the massage unit and start the unit.

3. Switch off the unit after use and store in a dry and safe place.


MAINTENANCE

Cleaning

1. Wipe the casing with a soft cloth, if necessary moistened with lukewarm water and a mild detergent. Do not use thinner, petrol, alcohol or the equivalent.
2. The electrodes can be cleaned with water if necessary. Allow to dry before use.
3. Store the unit dry, in a well-ventilated place.
4. Do not expose the unit to moisture, extreme temperatures or direct sunlight.

Jula reserves the right to make changes. In the event of problems, please contact our service department.

www.jula.com

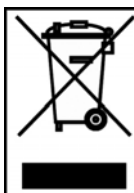Care for the environment!

Must not be discarded with household waste!
This product contains electrical or electronic components that should be recycled.
Leave the product for recycling at the designated station e.g. the local authority's recycling station.