# MARQUANT®

012808

## **SV** PERSONVÅG MED BLUETOOTH

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## SAFETY INSTRUCTIONS

- Do not use the product in damp or very hot areas.
- Temperature range: 5 to 35°C.
- Stand firmly on the product with your bare feet on the electrodes on the platform. Avoid putting weight near the edge of the scales to avoid the product tipping over, which could cause personal injury.
- Avoid dropping the product and do not subject it to knocks and blows – risk of material damage. This product is a precision instrument that must be handled carefully.
- Keep the product clean. Do not use strong or corrosive detergents.
- Do not use the product where there are strong vibrations. Stand firmly on the platform when in use.
- Measure your weight at the same time every day for best results.
- Replace the battery as soon as the battery charge is low, otherwise there is a risk that battery leakage will damage the product.
- Pregnant women, persons with oedema and persons who are undergoing dialysis should not use the product. Persons with pacemakers or other medical implants should not use the product.
- The calculated values for water content or body fat should only be used for reference and not for medical purposes. Consult a doctor if water content or body fat percentage is higher or lower than normal.
- Avoid contact with strong electromagnetic fields when using the product.
- Do not use the product on a wet surface.

## SYMBOLS

	Read the instructions.
	Approved in accordance with the relevant directives.
	Recycle discarded product in accordance with local regulations.

## TECHNICAL DATA

Max weight	180 kg
Size	L 30 x W 30 x H 2 cm

## DESCRIPTION

These scales with body fat measuring use a high precision strain gauge to register weight and changes in weight on your smartphone. The product can transfer weight data to a smartphone using Bluetooth, and create a curve of weight changes so that it is easy to see trends. The app must be downloaded in order to use the product.

## ASSEMBLY

### INSERTING THE BATTERIES

Open the battery compartment on the back and insert the batteries with the correct polarity.

Use 2x 1.5 V AAA batteries.

1. Remove the battery.
2. Put one end of the new battery under the battery holder and press down the battery in place. Make sure to put in the battery with the correct polarity.

**NOTE:**

**Replace the battery cover after inserting batteries. Only use batteries of the recommended type. Remove the batteries if the scales not going to be used for three months, or longer.**

**USE****DOWNLOADING AND INSTALLING APP**

Scan the QR code to download the Movinglife app. Or, search for "Movinglife" in the App Store or Google Play, and then download the app

**FIG. 1****Compatibility**

- Online version in App Store for IOS system.
- Android 4.3 or later.

**NOTE:**

**The app is continuously updated and therefore can differ from the information in these instructions. Search for updates to find the relevant information.**

**REGISTRATION AND CONNECTION**

Open the app MovingLife, connect the scales to the app and register your personal profile. Check before connecting that Bluetooth is enabled and that the display on the scales is on.

1. Step softly onto the scales to start them.
2. The display goes on.
3. Connect the scales to the app.

**FIG. 2****NOTE:**

**Make sure to specify your correct personal data to get correct results.**

**WEIGHING**

Place the scales on a hard, flat surface. Stand on the scales in your bare feet. The weight shown stabilises and then flashes, after which it is transferred. The scales switch off automatically after a few seconds. Data and history are stored and shown in the app.

1. Place the scales on a hard, flat surface.
2. Make sure that your bare feet are placed on the electrodes on the scales.
3. Check the data in the app.

**FIG. 3****NOTE:**

**To measure body fat correctly your bare feet must be placed on the electrodes on the scales.**

**BUTTON FUNCTIONS AND SETTING OF PERSONAL PARAMETERS**

SET	Start the scales and browse to next parameter.
▲	The value increases by one when pressed. Keep pressed to increase the values more quickly.
▼	The value decreases by one when pressed. Keep pressed to decrease the values more quickly.

**NOTE:**

1. **If the product model does not have any buttons, set personal parameters in the app.**
2. **The buttons are only intended for guests and the user code is as standard set to 0. It is not possible to set the user code on the scales. A user code is automatically allocated by the app. Set personal parameters in the app.**



The range of settings for the different parameters are given below.

<b>Gender</b>	<b>Height</b>	<b>Age</b>
Male/ female	100 to 250 cm	10 to 100 years

### **SETTING THE WEIGHT UNIT**

Set the required weight unit in the app or press the UNIT button on the back of the scales to switch weight unit.

## WARNINGS

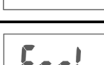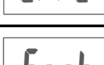Display	Significance	Explanation
	Body fat measuring.	Body fat measuring in progress. Stand still on the scales with your bare feet on the electrodes.
	Indication of low battery charge.	Low battery charging – replace the battery.
	Indication of overloading.	The weight on the scales exceeds the maximum permitted weight for the scales. Step off, otherwise the scales can be damaged.
	Calibration.	An error occurred during the weighing. Repeat the weighing.
	Low body fat.	Body fat percentage is at unhealthy low level.
	High body fat.	Body fat percentage is at unhealthy high level.

## TROUBLESHOOTING

Problem	Check/correct
The app does not receive the weighing result from the scales.	Check that the connection icon in the app shows that the smartphone is connected to the scales.
	Check that Bluetooth is enabled on the smartphone and that the scales and smartphone are within range.
	Check the batteries, and replace if necessary.
	Restart the app.
	Disable and enable Bluetooth, and repeat the connection procedure between the scales and the smartphone.
It is not possible to connect to the smartphone.	Check that only one smartphone is connected to the scales.
	Check that the smartphone has Bluetooth 4.0 or later.
Nothing is shown on the display.	Check the batteries, and replace if necessary.
	Repairs must only be carried out by an authorised service centre or qualified personnel.

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© Jula AB